



science
lab
research



innovation
quality
reputation
unsurpassed



Products ships in 5 days designated by a (*). Quantity restrictions apply—30 per table; 6 per cabinet.

ALPHABETICAL / CATEGORICAL

26	Accessories	40	Casework	39	Shower/Eyewash Station
16	ADA Workstations	23	Fume Hood Unit	13	Sink Units
10	Adjustable Height Tables	24	Instructor's Desks	42	Stools
36	Base Cabinets/Casework	33	Mirrors	7	Tables
9	Book Compartment Tables	28	Mobile Demonstration Labs	36	Tall Storage
41	Book Trucks/Bookcases	33	Mobile Storage	12	Wall Service Bench
10	Cabinet Tables	11	Riser Tables	14	Workstations

NUMERICAL ORDER

Page #	Item #	Product Name	Page #	Item #	Product Name
40	103-121 SERIES	BASE CABINETS	13	3116K	SERVICE ISLAND
24	1110 SERIES	8' INSTRUCTOR'S DESK	12	3200 SERIES	WALL SERVICE BENCH
24	1131K	TEACHER'S WORK DESK	13	3303K	CLEAN UP SINK
26	1140 SERIES	SIDE DESK	22	3400 & 3500 SERIES	PENINSULA WORKSTATIONS
25	1200 SERIES	5' INSTRUCTOR'S DESK	42	445-447 SERIES	CHEMICAL STORAGE UNITS
17	1500 & 1600 SERIES	OCTAGON WORKSTATIONS	33	4000/4001	MIRROR/MARKERBOARD
13	1536K	TRI-FACIAL SINK	33	4101K	MOBILE STORAGE CABINET
23	1800	FUME HOOD	32	4111 SERIES	MOBILE DEMO LAB
26	100002-7 SERIES	UPRIGHTS, CLAMPS, CROSSBARS	31	4121 SERIES	MOBILE DEMO LAB
27	100013-467 SERIES	FAUCETS, GAS COCKS	31	4222 SERIES	MOBILE DEMO LAB
8	100076X4	LEG BOOTS	30	4332 SERIES	MOBILE DEMO LAB
27	100118	WATER JUG	29	4342 SERIES	MOBILE DEMO LAB
36	100135	TOTE TRAYS	32	4352 SERIES	MOBILE DEMO LAB
26	100245GRARM	GREENLAW ARM	33	4400 SERIES	MOBILE STORAGE CABINET
40	100283	RUBBER BASE MOLDING	29	4500 SERIES	MOBILE DEMO CART
26	100604	UPRIGHT SET	28	4544	INTERACTIVE MEDIA CABINET
8	100728X4	FLOOR MOUNTING BRACKETS	35	4701K	MOBILE MICROSCOPE CABINET
26	100994	LABHAND	35	4741K	MOBILE MICRO-CHARGE STATION
9	2100-2200 SERIES	2-STUDENT CABINET TABLES	34	4751K	MOBILE TOTE STORAGE
21	2800 SERIES	PENINSULA WORKSTATIONS	34	4901K	WRITE N ROLL STORAGE
18	2894K	EZ-LIFT WORKSTATIONS	42	5000 SERIES	STOOLS
18	2900 SERIES	EZ-LIFT WORKSTATIONS	34	5201K	MOBILE BALANCE STORAGE
14	2916K, 2946K, 2956K	FORWARD VISION I	35	5301K	MOBILE LAPTOP STORAGE
16	2920 & 2960 SERIES	FORWARD VISION III	39	6900 SERIES	SHOWER/EYEWASH CABINETS
15	2944K	FORWARD VISION II	41	BT SERIES	BOOK TRUCKS
16	2946, 2946, ADA	FORWARD VISION ADA	9	C2100-C2200 SERIES	2-STUDENT CABINET TABLES
40	207000 SERIES	EPOXY RESIN TOPS	9	C2300-C2400 SERIES	4-STUDENT CABINET TABLES
40	214000 SERIES	PHENOLIC RESIN TOPS	20	C2400 SERIES	ISLAND WORKSTATIONS
23	235640	REPLACEMENT FUME HOOD FILTER	19	C2500 SERIES	ISLAND WORKSTATIONS
42	301 SERIES	BOOKCASE	20	C2600 SERIES	ISLAND WORKSTATIONS
36	313/315 SERIES	TALL STORAGE	19	C2700 SERIES	ISLAND WORKSTATIONS
36	350-4822	TOTE TRAY STORAGE CABINET	8	C7100-C7900	BOOK COMPARTMENT TABLES
37	351-4822	TOTE TRAY/SHELVING CABINET	40	D03-D06 SERIES	WALL CABINETS
36	353/358 SERIES	TALL STORAGE	8	L Series	REPLACEMENT TABLE PARTS
37	354-4830	ROCK/PAPER STORAGE CABINET	7	P7100-P7900	PLAIN APRON TABLES
38	360-3622	WARDROBE CABINET	10	P8000 SERIES	ADJUSTABLE HEIGHT TABLES
37	370-372 SERIES	MICROSCOPE STORAGE CABINET	11	PR7000 SERIES	TABLES WITH RISER
38	373-3616	MICRO-CHARGER CABINET	10	PR8000 SERIES	ADJ HT TABLE W/RISER
38	375-377 SERIES	SKELETON CABINETS			

www.diversifiedwoodcrafts.com

interactive



Image Downloads

Need an image? Multiple views? We have you covered...our online image library makes it easy to showcase our products.



Advanced Search

We have over 5,000 different tables. To find the product number that matches your needs, use our advanced table search. Find it super quick!



Product Cost

For dealers that are logged in to our site: do a product search; get all the product information; get the product cost. We make it simple!



Freight Cost

No need to call for a freight estimate. Use our system to get an up-to-date estimate 24/7.



Order Details

Find the order details...tracking number, when it was ordered, acknowledged and shipping information. Also find where it is going and any notes that would go on the "Bill of Lading."



Job Registration

If you have done the work and have us specified on a bid, just log-in to our web site and register that job.


Our Great Warranty

Products are guaranteed against material defect and/or faulty workmanship from plant for the lifetime of the product. During this period Diversified Woodcrafts will replace or repair, at its discretion, any product that upon inspection, exhibits defects in material or workmanship. This warranty shall include the costs of replacement parts or materials, but does not include labor, or shipping costs for such. This warranty doesn't cover damage in transit or damage caused by misuse, abuse or faulty installation, color grain or texture of wood or other covering materials. No other warranty is expressed or implied.



Our Stain & Finish Quality

All the products come in our Northwoods Oak stain, unless otherwise noted. Our finish is a chemical resistant, environmentally-friendly UV finish that meets strict SEFA standards, and is beautiful to boot!

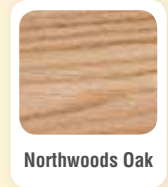

Our Customer Service

Our quality isn't just limited to our product. Dedicated customer service representatives can help you find the product that best fits your needs. Usually, you want it fast; so our top selling items are in our Fastrak program, and ship within five days of receipt of order. And you can count on our timely shipping...for five years in a row, we have had a 96% or better on-time shipping rate.


Our Shipping FAQs

- All our products are shipped: (FOB Wisconsin, 54174)
- Our shipping record...+96%!
- Table Shipping Cost: The weight of the pallet is included in the weight shown on the price list except in the case of tables. We ship up to ten tables per pallet, so adding the weight of a pallet to each individual table would cause the calculated weight to be excessive. Here is how it works:
 - Each square foot equals 5 lbs. of table top.
 - Take the number of tables and calculate how many pallets are needed
 - Example: A 24" x 60" table order of 18 tables would be 10 square feet, and equals 2 pallets. (10 sq. ft. x 5 lbs. = 50 lbs.); (50 lbs. x 2 pallets = 100 lbs)

our vitals



table guide

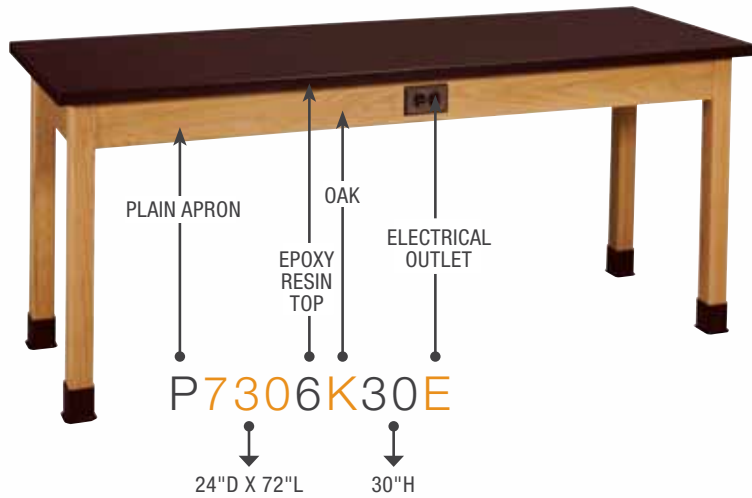

We offer a number of table options including tops, height and size. The coding system outlined below allows you to configure tables to meet your specific needs.

APRON STYLES	P = PLAIN	C = BOOK COMPARTMENT	D = DRAWER
Size	710 = 24" X 48" 711 = 42" X 48" 712 = 30" X 48" 713 = 30" X 54" 714 = 30" X 60" 715 = 30" X 72" 716 = 21" X 48"	717 = 36" X 48" 718 = 21" X 54" 719 = 36" X 54" 720 = 24" X 54" 721 = 21" X 60" 722 = 36" X 60" 723 = 21" X 72"	724 = 36" X 72" 730 = 24" X 72" 740 = 42" X 72" 748 = 48" Round 760 = 24" X 60" 780 = 42" X 54" 790 = 42" X 60"
Top	1 = 1-1/4" Plastic Laminate 2 = 1-1/4" ChemGuard 4 = 1" Phenolic Resin 5 = 1-3/4" Maple		
Wood	K = Oak M = Maple		
Height	30 = 30" H 34 = 34" H 36 = 36" H		
Options	N = None L = Leg Brace E = Electrical S = Leg Brace & Stretcher <small>(Recommended for tables 72" W. Required for tables over 72" W)</small>		
	LE = Leg Brace & Electrical SE = Leg Brace, Stretcher, Electrical C = Caster LC = Leg Brace, Caster SC = Leg Brace, Stretcher, Caster		



CORNER BRACES

Our table aprons have heavy-duty 13-gauge steel corner braces. They span the height of the apron and are dadoed and screwed into the apron. The leg is held in place with a hanger bolt, nut and lock washer to ensure a tight fit.



TOP ATTACHMENT

Tables are shipped with the tops attached by inserting screws through a minimum of ten pocket holes. This method is extremely strong and is essential for the table's stability.



LEG BOOTS

Our tables come with thick rubber leg boots that stay secure when the table is moved. These boots protect the legs from water and chemicals.



FINISH

Solid oak or maple (your choice) wood finished with a non-emitting UV finish that is extremely chemical resistant.

plain apron tables



PLAIN APRON TABLES—Are available in a variety of sizes, with several different top materials. These tables are built for years of dependable service. Legs and aprons are constructed of solid red oak, and are smoothly machined. They are finished with a chemical resistant, earth-friendly UV finish. Aprons are constructed of solid red oak and are 13/16"D x 4-1/2"H with heavy-duty 13-gauge steel corner braces. 2-1/4" square solid oak legs, with 5/16" hanger bolts installed, are attached to apron with nuts and washers supplied. All tables are equipped with non-skid adjustable glides for precise leveling and for prevention of rocking on uneven floors. Tables are supplied with black rubber leg boots which offer protection against scuffing and chemical spills. Standard table height is 30", with a weight capacity of 500 lbs. (Shipped KD with tops attached)



Size	P. Lam	ChemGuard	Solid Phenolic	Solid Epoxy	No Top
21" x 48"	P7161K30N	P7162K30N	P7164K30N	P7166K30N	P7160K30N
21" x 54"	P7181K30N	P7182K30N	P7184K30N	P7186K30N	P7180K30N
21" x 60"	P7211K30N	P7212K30N	P7214K30N	P7216K30N	P7210K30N
21" x 72"	P7231K30N	P7232K30N	P7234K30N	P7236K30N	P7230K30N
24" x 48"	P7101K30N*	P7102K30N*	P7104K30N	P7106K30N*	P7100K30N
24" x 54"	P7201K30N*	P7202K30N*	P7204K30N	P7206K30N*	P7200K30N
24" x 60"	P7601K30N*	P7602K30N*	P7604K30N	P7606K30N*	P7600K30N
24" x 72"	P7301K30N*	P7302K30N*	P7304K30N	P7306K30N*	P7300K30N
30" x 48"	P7121K30N	P7122K30N	P7124K30N	P7126K30N	P7120K30N
30" x 54"	P7131K30N	P7132K30N	P7134K30N	P7136K30N	P7130K30N
30" x 60"	P7141K30N*	P7142K30N*	P7144K30N	P7146K30N*	P7140K30N
30" x 72"	P7151K30N	P7152K30N*	P7154K30N	P7156K30N	P7150K30N
36" x 48"	P7171K30N	P7172K30N	P7174K30N	P7176K30N	P7170K30N
36" x 54"	P7191K30N	P7192K30N	P7194K30N	P7196K30N	P7190K30N
36" x 60"	P7221K30N	P7222K30N	P7224K30N	P7226K30N	P7220K30N
36" x 72"	P7241K30N	P7242K30N	P7244K30N	P7246K30N	P7240K30N
42" x 48"	P7111K30N	P7112K30N	P7114K30N	P7116K30N	P7110K30N
42" x 54"	P7801K30N	P7802K30N	P7804K30N	P7806K30N	P7800K30N
42" x 60"	P7901K30N	P7902K30N	P7904K30N	P7906K30N	P7900K30N
42" x 72"	P7401K30N	P7402K30N	P7404K30N	P7406K30N	P7400K30N
48" Dia.	P7481K30N	P7482K30N	P7484K30N	P7486K30N	P7480K30N



book compartment tables



BOOK COMPARTMENT TABLES—Are built with the same quality techniques as the plain apron table. They are finished with a chemical resistant, earth-friendly UV finish. These tables have either two or four book compartments (depending on size) installed in the apron. The book compartments measure 14-1/4"W x 12-3/4"D x 3-7/8"H, with a weight capacity of 500 lbs. On tables 72"W or over, we recommend ordering with leg braces and stretchers. (Shipped KD with tops attached)

Size	Compartments	P. Lam	ChemGuard	Solid Phenolic	Solid Epoxy	No Top
21" x 48"	2	C7161K30N	C7162K30N	C7164K30N	C7166K30N	C7160K30N
21" x 54"	2	C7181K30N	C7182K30N	C7184K30N	C7186K30N	C7180K30N
21" x 60"	2	C7211K30N	C7212K30N	C7214K30N	C7216K30N	C7210K30N
21" x 72"	2	C7231K30N	C7232K30N	C7234K30N	C7236K30N	C7230K30N
24" x 48"	2	C7101K30N*	C7102K30N*	C7104K30N	C7106K30N*	C7100K30N
24" x 54"	2	C7201K30N*	C7202K30N*	C7204K30N	C7206K30N*	C7200K30N
24" x 60"	2	C7601K30N*	C7602K30N*	C7604K30N	C7606K30N*	C7600K30N
24" x 72"	2	C7301K30N*	C7302K30N*	C7304K30N	C7306K30N*	C7300K30N
30" x 48"	2	C7121K30N	C7122K30N	C7124K30N	C7126K30N	C7120K30N
30" x 54"	2	C7131K30N	C7132K30N	C7134K30N	C7136K30N	C7130K30N
30" x 60"	2	C7141K30N*	C7142K30N*	C7144K30N	C7146K30N*	C7140K30N
30" x 72"	2	C7151K30N	C7152K30N*	C7154K30N	C7156K30N	C7150K30N
36" x 48"	2	C7171K30N	C7172K30N	C7174K30N	C7176K30N	C7170K30N
36" x 54"	2	C7191K30N	C7192K30N	C7194K30N	C7196K30N	C7190K30N
36" x 60"	2	C7221K30N	C7222K30N	C7224K30N	C7226K30N	C7220K30N
36" x 72"	2	C7241K30N	C7242K30N	C7244K30N	C7246K30N	C7240K30N
42" x 48"	4	C7111K30N	C7112K30N	C7114K30N	C7116K30N	C7110K30N
42" x 54"	4	C7801K30N	C7802K30N	C7804K30N	C7806K30N	C7800K30N
42" x 60"	4	C7901K30N	C7902K30N	C7904K30N	C7906K30N	C7900K30N
42" x 72"	4	C7401K30N	C7402K30N	C7404K30N	C7406K30N	C7400K30N

REPLACEMENT PARTS

* QUANTITY RESTRICTIONS APPLY

Unit #	Description	Unit #	Description
L30-REP*	30"H Leg (1 each)	L3036*	36"H Leg (set of 4)
L34-REP	34"H Leg (1 each)	100076X4*	Leg Boots (set of 4)
L36-REP*	36"H Leg (1 each)	100728X4*	Floor Mounting Brackets
L3030*	30"H Leg (set of 4)	(set of 4)	
L3034	34"H Leg (set of 4)		

cabinet tables



STORAGE CABINET TABLES—Are constructed of solid oak and hardwood veneers, this table features a center cabinet with lock, and rubber base molding. This table is available with plain or book compartment aprons and comes in both single (2-student) and double faced (4-student) models. The center cabinet is available with two locking drawers or as a plain base cabinet with adjustable shelf. They are finished with a chemical resistant, earth-friendly UV finish. This unit is available with our complete line of top materials to satisfy every requirement or budget.



Single Face (2-Student)	Wt.	Double Face (4-Student)	Wt.	Apron	Cabinet	Top
2101K	265 lbs.	2301K	380 lbs.	Plain	Plain	P. Lam
C2101K	275 lbs.	C2301K	400 lbs.	Compartment	Plain	P. Lam
2201K	275 lbs.	2401K	400 lbs.	Plain	Drawers	P. Lam
C2201K	275 lbs.	C2401K	410 lbs.	Compartment	Drawers	P. Lam
2102K	265 lbs.	2302K	380 lbs.	Plain	Plain	ChemGuard
C2102K	275 lbs.	C2302K	400 lbs.	Compartment	Plain	ChemGuard
2202K	275 lbs.	2402K	400 lbs.	Plain	Drawers	ChemGuard
C2202K	275 lbs.	C2402K	410 lbs.	Compartment	Drawers	ChemGuard
2104K	316 lbs.	2304K	443 lbs.	Plain	Plain	Solid Phenolic
C2104K	326 lbs.	C2304K	463 lbs.	Compartment	Plain	Solid Phenolic
2204K	330 lbs.	2404K	463 lbs.	Plain	Drawers	Solid Phenolic
C2204K	335 lbs.	C2404K	468 lbs.	Compartment	Drawers	Solid Phenolic
2106K	340 lbs.	2306K	500 lbs.	Plain	Plain	Solid Epoxy
C2106K	350 lbs.	C2306K	520 lbs.	Compartment	Plain	Solid Epoxy
2206K	350 lbs.	2406K	520 lbs.	Plain	Drawer	Solid Epoxy
C2206K	355 lbs.	C2406K	525 lbs.	Compartment	Drawers	Solid Epoxy



SINGLE FACE UNIT DIMENSIONS: (66"W x 24"D x 30"H)

- 1 Single-Faced Cabinet: 24"W x 22"D x 29"H
- 2 Book Compartments: 12"W x 13"D x 3"H
- 2 Drawers: 11-1/2"W x 18"D x 5-1/2"H
- Cabinet Bases w/full door have an adjustable shelf, with a 40 lb. capacity.



DOUBLE FACE UNIT DIMENSIONS: (66"W x 42"D x 30"H)

- 1 Double-Faced Cabinet: 24"W x 40"D x 29"H
- 4 Book Compartments: 12"W x 13"D x 3"H
- 4 Drawers: 11-1/2"W x 18"D x 5-1/2"H
- Cabinet Bases w/full door have an adjustable shelf, with a 40 lb. capacity.



adjustable height tables



ADJUSTABLE HEIGHT TABLE—A perfect solution for ADA requirements with the ability to adjust in 2 cm increments from 27"H to 39"H. Can also be adjusted for use as an incline surface for physics. Solid 4-1/2" oak aprons are attached to your choice of three top options. The wood is finished with a chemical resistant, earth-friendly UV finish. Optional casters sold separately.

Size	P. Lam	ChemGuard	Solid Phenolic
24" x 48"	P8101K*	P8102K*	P8104K
24" x 54"	P8201K*	P8202K*	P8204K
24" x 60"	P8601K*	P8602K*	P8604K
24" x 72"	P8301K*	P8302K*	P8304K
30" x 60"	P8141K*	P8142K*	P8144K
Casters Option	100709X4*	(Set of 4)	

* QUANTITY RESTRICTIONS APPLY



QUICK HEIGHT ADJUSTMENT
Adjusted in less than 30 seconds.



CASTERS
Quick release casters. Easy to attach.

new!

ADJUSTABLE HEIGHT RISER TABLE—This table has everything you could want. It adjusts easily from 27" to 39"H in 2 cm increments in less than 30 seconds without the use of tools. Its apron is solid oak and finished with a chemical resistant finish. It has a riser that provides extra workspace for each person to put his or her tools or laptop. The table has three top options: plastic laminate, chemical resistant laminate and phenolic resin. The riser's horizontal surface matches your choice of top surface. The 12"D riser is supported by six 1" diameter powder coated steel tubes. The riser is 10" above the table surface. Tabletop measures 30"D x 60"W.

Item #	Top
PR8141K	P. Lam
PR8142K	Chemical Resistant Laminate
PR8144K	Phenolic Resin



riser tables



PR7144K30N (Oak)

new!

FIXED HEIGHT RISER TABLE—Perfect table for STEM projects. We've taken our popular oak and maple science tables and enhanced them by adding a riser. The table features solid oak or maple (your choice) aprons and riser supports finished with a chemical resistant finish. The riser provides 10" of clearance above the work surface, is 12"D and made of the same material as the top you choose. Top options include: plastic laminate, chemical resistant laminate, phenolic resin and epoxy resin. 30"D and 36"D tables have the riser positioned at the back of the table. On 42"D tables the riser is centered on the table. All 72"W tables include leg braces and a stretcher for added strength. Tables stand 30"H. The top of the riser is 41" from the floor.



PR7144M30N (Maple)



PR7144K30N (Oak Back)

Size	P. Lam Oak	P. Lam Maple	Chem Oak	Chem Maple	Phenolic Oak	Phenolic Maple	Epoxy Oak	Epoxy Maple
30x48	PR7121K30N	PR7121M30N	PR7122K30N	PR7122M30N	PR7124K30N	PR7124M30N	PR7126K30N	PR7126M30N
30x54	PR7131K30N	PR7131M30N	PR7132K30N	PR7132M30N	PR7134K30N	PR7134M30N	PR7136K30N	PR7136M30N
30x60	PR7141K30N	PR7141M30N	PR7142K30N	PR7142M30N	PR7144K30N	PR7144M30N	PR7146K30N	PR7146M30N
30x72	PR7151K30S	PR7151M30S	PR7152K30S	PR7152M30S	PR7154K30S	PR7154M30S	PR7156K30S	PR7156M30S
36x48	PR7171K30N	PR7171M30N	PR7172K30N	PR7172M30N	PR7174K30N	PR7174M30N	PR7176K30N	PR7176M30N
36x54	PR7191K30N	PR7191M30N	PR7192K30N	PR7192M30N	PR7194K30N	PR7194M30N	PR7196K30N	PR7196M30N
36x60	PR7221K30N	PR7221M30N	PR7222K30N	PR7222M30N	PR7224K30N	PR7224M30N	PR7226K30N	PR7226M30N
36x72	PR7241K30S	PR7241M30S	PR7242K30S	PR7242M30S	PR7244K30S	PR7244M30S	PR7246K30S	PR7246M30S
42x54	PR7801K30N	PR7801M30N	PR7802K30N	PR7802M30N	PR7804K30N	PR7804M30N	PR7806K30N	PR7806M30N
42x60	PR7901K30N	PR7901M30N	PR7902K30N	PR7902M30N	PR7904K30N	PR7904M30N	PR7906K30N	PR7906M30N
42x72	PR7401K30S	PR7401M30S	PR7402K30S	PR7402M30S	PR7404K30S	PR7404M30S	PR7406K30S	PR7406M30S



wall service benches



3216K (Shown w/two 36"H Plain Apron Tables)



3216K



3226K



3236K



3246K

WALL SERVICE BENCH—Need all your services in one place or looking to add a wall of cabinets that include services? The wall service bench offers storage and services all in one unit. Constructed of solid hardwood and oak veneer plywood this unit features your choice of a 1" solid epoxy resin or phenolic resin top with a 18"W x 15"D x 11"H epoxy resin sink and two combination water/gas fixtures. The sink area allows access to the plumbing and provides a large storage space with two GFI protected AC duplex electrical outlets. There is enough room in the sink area to allow two 24"D tables to be butted against the unit for use as a peninsula work area. The outer ends of the unit offer four storage solutions:

- Doors with Adjustable Shelves
- Doors/Drawer Combination
- 4 Drawers/Doors
- 4 Drawers

All drawers and doors have locks and ADA compliant pulls. The unit is finished with a chemical resistant, earth-friendly UV finish and includes rubber base molding. Dimensions: 108"L x 24"D x 36"H.

Epoxy	Phenolic	Cabinet Style
3216K	3214K	Door
3226K	3224K	Door/Drawer
3236K	3234K	Door/4 Drawer
3246K	3244K	Drawer

sink units



TRIFACIAL WORKSTATION—Hexagonal workstation designed to be used with standard laboratory tables. This unit features an epoxy resin top with a drop-in hexagonal lipped sink. Apron and pedestal are constructed of solid hardwoods, and oak veneers, with a chemical resistant, earth-friendly UV finish. Equipped with three combination water/gas fixtures and three GFI duplex receptacles. The center pedestal has a removable panel to access utilities. Unit is also available in a flat top version without the sink and fixtures.

Item #	Size	Top
1536K	49"W x 56-1/2"D x 36"H	Solid Epoxy
1536KF	49"W x 56-1/2"D x 36"H w/Flat Epoxy Top	Solid Epoxy



SERVICE ISLAND—A compact service island constructed of solid oak, hardwood and oak veneers, with a chemical resistant, earth-friendly UV finish. This unit can be grouped with tables to provide an efficient working arrangement. The unit includes one 14" x 10" x 6" solid epoxy resin sink, two combination water/gas fixtures, two GFI AC duplex electrical receptacles, chemical-resistant polypropylene trap, and locking cupboard.

Item #	Size	Top
3116K	26"W x 24"D x 36"H	Solid Epoxy

CLEAN UP SINK—One piece top including sink is molded of polyolefin. (Polyolefin is easy to clean, extremely chemical resistant, and is highly impact resistant to reduce glassware breakage). Long grooved sloping counter leads to large sink in center. The low temperature strength eliminates cracking with dry ice, and the sink measures (24-1/2"W x 18-1/2"D x 14"H). It includes a hot and cold water mixing faucet. The cabinet base is constructed of hardwood and oak veneers with a chemical resistant, earth-friendly UV finish. Two adjustable shelves can be found inside the outer cabinets. It includes pulls, locks and rubber base molding.

Item #	Size	Top
3303K*	55-1/2"W x 28"D x 36-1/2"H	Polyolefin





workstations



2946K | Front View



2936K | Front View



2936KLH | Front View



2946K | Back View

FORWARD VISION I—This popular workstation meets the needs of students and teachers today. The Forward Vision 4-Student Workstation allows four students to work at a single station while all are facing the teacher. The unique design leaves enough room to store four stools (not included) under the workstation. Sink units comes with two combination water/gas fixtures. The 2916KLH and 2936KLH include two rod sockets and LabHands for holding laptop computers and books. All the units come with two GFI electrical outlets, and ample storage. Choose from two models of cabinets: Door cabinet with two adjustable shelves or the door/drawer cabinet with one adjustable shelf and a drawer. The stations are finished with a chemical resistant, earth-friendly UV finish. All doors and drawers lock. Dimensions: 96"W x 50"D x 36"H.

Style	Sink Dimensions	Door Cabinet	Door/Drawer Cabinet	Top
End Sink	22-3/8" x 15" x 17"	2916K	2936K	Solid Epoxy
End Sink w/LabHands	22-3/8" x 15" x 17"	2916KLH	2936KLH	Solid Epoxy
Drop-In Sink	12" x 8" x 6"	2946K*	2956K	Solid Epoxy
No Sink or Fixtures	N/A	2946KF	2956KF	Solid Epoxy

workstations



2944K | Front View



2944K | Back View

FORWARD VISION II—Everything science teachers want. Designed with teacher input and desire for a safe, secure classroom. Teachers maintain control of the class with all students facing the instructor. The large work surface—46"D x 88"W—is made of chemical resistant phenolic resin and holds the solid epoxy sink, and two combination water & gas fixtures. The cabinets provide locked storage space. Side mounted hooks hold book bags off the floor. The units are finished with a chemical resistant, earth-friendly UV finish. Includes two GFI AC duplex units. Also available with a flat top (no sink or fixtures). This model is compact enough to fit multiple units in a small classroom and is easy to assemble. Dimensions 88"W x 46"D x 36"H. Wt. 503 lbs.

Style	Item #	Sink Dimensions	Top
Forward Vision II	2944K*	12" x 8" x 6"	Solid Phenolic
Forward Vision II w/Flat Top	2944KF	N/A	Solid Phenolic



workstations



FORWARD VISION III—Durable lab table features a large 68"x 36" tabletop. It comes standard with two combination water/gas fixtures and one GFI AC electrical receptacle. Your choice of two cabinets: Full door cabinet with adjustable shelves or door/drawer cabinet with a drawer and storage area with one adjustable shelf. The door/drawer cabinet is equipped with a drawer, and a cabinet with one adjustable shelf. They are finished with a chemical resistant, earth-friendly UV finish. An 18" wide removable panel is located between the storage cabinets and provides access to all utilities. Unit can be ordered with a drop-in sink or with a flat top without sink or fixtures. Assembly and installation of electricity and plumbing is required. Dimensions: 68"W x 36"D x 36"H.

Sink	Sink Dimensions	Door Cabinet	Door/Drawer Cabinet	Top
Yes	12" x 8" x 6"	2926K	2966K	Solid Epoxy
No	N/A	2926KF	2966KF	Solid Epoxy
Yes	12" x 8" x 6"	2924K	2964K	Solid Phenolic
No	N/A	2924KF	2964KF	Solid Phenolic

new!

FORWARD VISION ADA —We've taken our most popular workstation and have made it ADA compatible. Constructed of hardwoods and oak veneers and finished with a chemical resistant finish. The 1" epoxy resin top is available with or without a sink and two combination water/gas fixtures with paddle handles. The sink includes the trap, strainer and stopper. The sloping front allows a person in a wheelchair to roll right up to the sink area. The door/drawer cabinet features a dovetailed drawer and storage area with an adjustable shelf, locks, and ADA compliant pulls. The door cabinet option includes everything previously listed except the drawer and includes 2 adjustable shelves. The unit also includes 2 GFI protected duplex outlets and rubber base molding. Unit measures: 96"W x 50"D x 34"H.



Item #	Cabinet	Sink Dimensions
2946K-L-ADA	Workstation with door cabinet	15"W x 18"L x 5"D
2946KF-L-ADA	Workstation with door cabinet w/flat top	N/A
2956K-L-ADA	Workstation with door/drawer cabinet	15"W x 18"L x 5"D
2956KF-L-ADA	Workstation with door/drawer cabinet w/flat top	N/A

workstations



Drawer Base Option



Pedestal Base Option

OCTAGON WORKSTATION—Crafted of solid hardwood/oak veneers, and finished with a chemical resistant, earth-friendly UV finish. Both units provide ample space for 4 people and are available in 56" or 62", solid 1" epoxy resin or phenolic resin tops. The pedestal with apron comes equipped with 4 drawers, the drawer cabinet unit comes equipped with 16 drawers. Both units come with 4 GFI AC duplex receptacles, and can be ordered with these options: two combination water/gas fixtures, 4 rod sockets and, a 16" x 16" x 7" solid epoxy sink with strainer, stopper and trap. Easily accessible to the utilities. All units are 36"H. (Shipped KD; FOB Wisconsin, 54174) Drawer Inside Dimensions:

- 1516K includes 4 each (18"W x 11"D x 3-1/2"H)
- 1546K includes 16 each (20-1/2"W x 11"D x 5-1/2"H)
- 1616K includes 4 each (24"W x 11"D x 3-1/2"H)
- 1646K includes 16 each 26-1/2"W x 11"D x 5-1/2"H

Style	Size	Top	With Sink	Flat Top*
Pedestal	56"W	Solid Epoxy	1516K	1516KF
Drawer	56"W	Solid Epoxy	1546K	1546KF
Pedestal	62"W	Solid Epoxy	1616K	1616KF
Drawer	62"W	Solid Epoxy	1646K	1646KF
Pedestal	56"W	Solid Phenolic	1514K	1514KF
Drawer	56"W	Solid Phenolic	1544K	1544KF
Pedestal	62"W	Solid Phenolic	1614K	1614KF
Drawer	62"W	Solid Phenolic	1644K	1644KF



workstations



EZ-LIFT WORKSTATIONS (MANUAL & PUSH-BUTTON)—

Both of the units' height is adjustable from 28" to 40"H. The students have plenty of space with either the 96"W x 50"D or 88"W x 46" work surface. Two hooks mounted underneath hold up to four book bags. The bone powder-coated steel and oak veneer base conceals the lifting mechanism. It is finished with a chemical resistant, earth-friendly UV finish. The unit provides six GFI protected electrical outlets. Assembly and power source required.

Manual-Crank: Fold-away crank accessible on the student side of the workstation.

Push-Button: Easy, power option. The motor is accessible by removing one of the oak panels. Plug only supplies power to the hydraulics on the push-button model. AC Duplex needs to be hard-wired in the field.



Fold-Away Crank


Item #	Top Size	Top	Adjustable Unit
2994K	96"W x 50"D	Solid Phenolic	Manual-Crank
2996K	96"W x 50"D	Solid Epoxy	Manual-Crank
2904K	88"W x 46"D	Solid Phenolic	Manual-Crank
2906K	88"W x 46"D	Solid Epoxy	Manual-Crank
2974K	96"W x 50"D	Solid Phenolic	Push-Button
2976K	96"W x 50"D	Solid Epoxy	Push-Button
2984K	88"W x 46"D	Solid Phenolic	Push-Button
2986K	88"W x 46"D	Solid Epoxy	Push-Button

EZ-LIFT WORKSTATIONS W/STORAGE—Our most flexible workstation yet. This unit is designed to be used by four students, all facing the instructor. The fold-away crank allows the units height to be adjusted from 28"H to 40"H, for ease of use with a wheelchair. The front of the unit has locking storage measuring 36"W x 8"D x 24"H. Inside the cabinet are two lipped shelves and a quad electrical outlet. This space can be used for computer storage, recharging or for general storage. The front panel includes two more duplex outlets and two data ports. The unit is finished with a chemical resistant, earth-friendly UV finish, and boosts 500 lbs. weight capacity. The top is 1" solid phenolic resin. Overall size: 96"W x 54"D x 28"-40"H. Assembly and power source required.



Item #	Size	Top	Adjustable Unit
2894K	96"W x 54"D x 28"- 40"H	Solid Phenolic	Manual-Crank

workstations

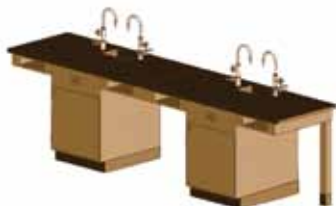

Two-Student Service Centers are ganged together to form a continuous row of workstations. This example illustrates the appearance of the C2716K. Two-Station Service Center when it is expanded to become the C2726K 4-Student Workstation.



TWO-TO-EIGHT STUDENT SERVICE ISLAND—Single faced tables that provide a fully equipped science center from two to eight students. Each unit is constructed of solid oak, and hardwood veneers. The single faced units contain a solid epoxy resin sink (12" x 8" x 6"H) along with faucets for every two students. Every unit contains enough gas cocks for each student along with GFI AC duplex electrical receptacles, rod sockets and book compartments. All units contain locks for the drawers and doors. Choose from models with a full cupboard only or a half cupboard and drawers. Can be ordered without sink, fixtures and rod sockets. They are finished with a chemical resistant, earth-friendly UV finish. Single face two-student unit measures 66"L x 30"W x 36"H. Add 66" in length for each additional two students. Your choice of 1" solid epoxy resin or 1" phenolic resin tops.

w/Sink Fixtures & Rod Sockets	w/o Sink Fixtures & Rod Sockets	# Of Students	Top	Cabinet Style
C2514K	C2514KF	2	Solid Phenolic	Cupboard w/ Drawers
C2524K	C2524KF	4	Solid Phenolic	Cupboard w/ Drawers
C2534K	C2534KF	6	Solid Phenolic	Cupboard w/ Drawers
C2544K	C2544KF	8	Solid Phenolic	Cupboard w/ Drawers
C2516K	C2516KF	2	Solid Epoxy	Cupboard w/ Drawers
C2526K	C2526KF	4	Solid Epoxy	Cupboard w/ Drawers
C2536K	C2536KF	6	Solid Epoxy	Cupboard w/ Drawers
C2546K	C2546KF	8	Solid Epoxy	Cupboard w/ Drawers
C2714K	C2714KF	2	Solid Phenolic	Cupboard
C2724K	C2724KF	4	Solid Phenolic	Cupboard
C2734K	C2734KF	6	Solid Phenolic	Cupboard
C2744K	C2744KF	8	Solid Phenolic	Cupboard
C2716K	C2716KF	2	Solid Epoxy	Cupboard
C2726K	C2726KF	4	Solid Epoxy	Cupboard
C2736K	C2736KF	6	Solid Epoxy	Cupboard
C2746K	C2746KF	8	Solid Epoxy	Cupboard



workstations



C2416K



Four-Student Service Centers are ganged together to form a continuous row of workstations. This example illustrates the appearance of the C2416K Four-Station Service Center when it is expanded to the 2426K Eight-Student Workstation.

FOUR-TO-SIXTEEN STUDENT

SERVICE ISLAND—Double faced tables that provide a fully equipped science center from four to sixteen students. Each unit is constructed of solid oak, oak and hardwood veneers and finished with a chemical resistant, earth-friendly UV finish. The double faced units contain a solid epoxy resin sink (12" x 8" x 6"H) along with faucets for every two students. Every unit contains enough gas cocks for each student along with GFI AC duplex electrical receptacles, rod sockets and book compartments. All units contain locks for the drawers and doors. Choose from models with a full cupboard only or a half cupboard and drawers. Can be ordered without sink, fixtures and rod sockets. Double face four student unit measures 66"L x 48"W x 36"H. Add 66" in length for each additional four students.

w/Sink Fixtures & Rod Sockets	w/o Sink Fixtures & Rod Sockets	# Of Students	Top	Cabinet Style
C2414K	C2414KF	4	Solid Phenolic	Cupboard w/Drawers
C2424K	C2424KF	8	Solid Phenolic	Cupboard w/Drawers
C2434K	C2434KF	12	Solid Phenolic	Cupboard w/Drawers
C2444K	C2444KF	16	Solid Phenolic	Cupboard w/Drawers
C2416K	C2416KF	4	Solid Epoxy	Cupboard w/Drawers
C2426K	C2426KF	8	Solid Epoxy	Cupboard w/Drawers
C2436K	C2436KF	12	Solid Epoxy	Cupboard w/Drawers
C2446K	C2446KF	16	Solid Epoxy	Cupboard w/Drawers
C2614K	C2614KF	4	Solid Phenolic	Cupboard
C2624K	C2624KF	8	Solid Phenolic	Cupboard
C2634K	C2634KF	12	Solid Phenolic	Cupboard
C2644K	C2644KF	16	Solid Phenolic	Cupboard
C2616K	C2616KF	4	Solid Epoxy	Cupboard
C2626K	C2626KF	8	Solid Epoxy	Cupboard
C2636K	C2636KF	12	Solid Epoxy	Cupboard
C2646K	C2646KF	16	Solid Epoxy	Cupboard

workstations



PERIMETER WORKSTATIONS—Provides work space for up to four-students while storing your laboratory equipment in these perimeter workstations. Choose from 4 different configurations (see below). The four cabinets which make up a single unit measure 24"W x 16"D x 35"H. The stations are designed to utilize space along walls where they can stand alone or be ganged together. Durable solid oak and oak veneer construction finished with a chemical resistant, earth-friendly UV finish. Each unit includes: your choice of 1" solid epoxy resin or 1" solid phenolic resin top with back splash, solid epoxy drop-in style sink with strainer, 2 combination water and gas fixtures, 4 rod sockets, rod storage (rods not included), 2 GFI electrical receptacles, keyed-alike locks on all doors and drawers and rubber base molding. Sink measures: 12"W x 8"D x 6"H. Total unit measure 90"W x 84"D x 36"H. Work surface on wall portion is 18"D.



2826K | Front View



2816K | Front View



2836K | Front View



2846K | Front View

Item #	w/Sink & Fixtures	Top	Cabinet Style
2814K	Yes	Solid Phenolic	Door
2824K	Yes	Solid Phenolic	Door/Drawer Unit
2834K	Yes	Solid Phenolic	Door/4 Drawer Unit
2844K	Yes	Solid Phenolic	4 Drawer Unit
2816K	Yes	Solid Epoxy	Door
2826K	Yes	Solid Epoxy	Door/Drawer Unit
2836K	Yes	Solid Epoxy	Door/4 Drawer Unit
2846K	Yes	Solid Epoxy	4 Drawer Unit
2814KF	No	Solid Phenolic	Door
2824KF	No	Solid Phenolic	Door/Drawer Unit
2834KF	No	Solid Phenolic	Door/4 Drawer Unit
2844KF	No	Solid Phenolic	4 Drawer Unit
2816KF	No	Solid Epoxy	Door
2826KF	No	Solid Epoxy	Door/Drawer Unit
2836KF	No	Solid Epoxy	Door/4 Drawer Unit
2846KF	No	Solid Epoxy	4 Drawer Unit
100604	Optional Upright Rod Set w/Full Crossbar		



workstations




3416K-R (Sink in Peninsula)



3516K-R (Sink in Wall)



3516KF-R (No Sink)

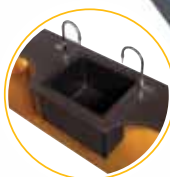

new!

PENINSULA WORKSTATION—This unique L-shaped unit is designed so four students can use the workstation while facing the teacher. The unit has storage on the peninsula and on the wall. It is made of oak and oak veneers with black powder-coated metal accents. The unit can be ordered with either a 1" powdercoated or epoxy resin top. Three sink options are available: epoxy sink (18"L x 15"W x 7.9"D) on the wall section, epoxy sink (16"L x 8"W x 6.8"D) in the peninsula or a no sink option. Units with sinks include cold water fixtures (call for quote if gas is needed). All units include a pop-up duplex electrical outlet. Finished in natural oak with an earth-friendly chemical resistant UV finish. Wall portion measures: 96"W x 12"D; peninsula extends: 96"; overall size: 110"D x 96"W x 36"H.

Item #	Sink Location	Top	Room Orientation
3416K-L	Sink in Peninsula	Epoxy	Left Side
3416K-R	Sink in Peninsula	Epoxy	Right Side
3419K-L	Sink in Peninsula	Powdercoat	Left Side
3419K-R	Sink in Peninsula	Powdercoat	Right Side
3516K-L	Sink in Wall	Epoxy	Left Side
3516K-R	Sink in Wall	Epoxy	Right Side
3519K-L	Sink in Wall	Powdercoat	Left Side
3519K-R	Sink in Wall	Powdercoat	Right Side
3516KF-L	No Sink	Epox	Left Side
3516KF-R	No Sink	Epoxy	Right Side
3519KF-L	No Sink	Powdercoat	Left Side
3519KF-R	No Sink	Powdercoat	Right Side

Room orientation is from the perspective of the student facing the teacher.

fume hood unit



1800K (Back)

new!

MOBILE FUME HOOD STATION—A versatile fume hood that can be used in multiple classrooms. This ductless fume hood has an all-inclusive filter that handles liquids, gasses and solid chemicals. This unit is so flexible, it does not require ducting (saves energy) and plugs into an ordinary outlet. The fume hood part of the unit is enclosed by 5/16" thick acrylic so a teacher can demonstrate while students surround the hood. The metallic parts are anti-corrosion metallic alloy protected by a thermo-hardened anti-acid polymer coating. The internal hood dimensions are 30-1/8"W x 21-1/4"D x 33-3/4"H; external dimensions are 31-1/2"W x 25"D x 45"H. Overall dimensions are: 33-3/8"W x 28"D x 81"H. Unit comes with an air flow meter, flow monitor, sampling port, programmable timer, energy ports, one filter (2 year average life) and list of approved chemicals. Fume hood comes standard with an internal light fixture. The fume hood is mounted to a 28"D x 32-3/8"W x 36"H oak cabinet. It provides an easy way to move the fume hood within a room or school. It rolls easily on 4" casters and has two drawers, two pull-out lipped shelves lined with a chemical resistant liner. Unit is finished in natural oak with an earth-friendly chemical resistant UV finish.



Pull-out Shelves




Chemical Monitor



Filter

Item #	Product
1800K	Mobile Fume Hood Station (Includes Filter and Cabinet)
235640	Replacement filter

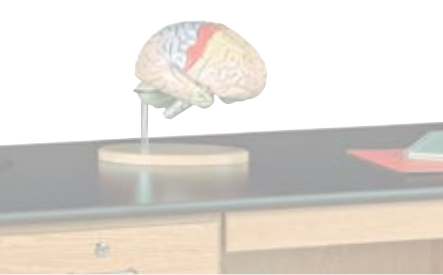

instructor's desks



Sink can be located on the left or right.

8' INSTRUCTOR'S DESK—Constructed of solid oak, oak and hardwood veneers, this unit serves as a complete demonstration center offering a spacious work surface. It is finished with a chemical resistant, earth-friendly UV finish. Included are drawers, and a cupboard for storage of large items. The drawers are expertly crafted from the finest materials with dovetail joinery throughout. Drawers run on epoxy-coated steel glides with nylon bearings for smooth operation. Equipment includes: one mixing faucet, double gas cock, GFI protected AC duplex receptacle, polypropylene trap, set of support rods and crossbar, rubber base molding and locks for cupboard and all five drawers (two small, three large). Can also be ordered without the sink, faucet, gas cock, rod sockets or support rod. Note: mirror, and mirror/markerboard will not fit this unit. Dimensions:

- Knee space: 22"W x 32-1/2"H,
- Overall: 96"L x 30"D x 36"H.
- Small Drawers: 14-3/4"W x 19"D x 5-1/2"H
- Large Drawers: 33-3/4"W x 19"D x 5-1/2"H
- Cupboard: 34-1/2"W x 28"D x 29-1/4"H
- Sink: 18"W x 15"D x 11"H

Item #	With Sink	Top
1114K	Yes	Solid Phenolic
1116K*	Yes	Solid Epoxy
1114KF	No	Solid Phenolic
1116KF	No	Solid Epoxy

TEACHER'S WORK DESK—A beautiful desk that will match the rest of your furniture. This 60"W x 30"D x 30"H desk is constructed of select hardwoods and oak veneers. It is finished with a chemical resistant, earth-friendly UV finish. The top is 1" black high pressure laminate. The top also has two grommet holes for threading power cords. The pedestals have 2 file drawers, four smaller drawers and one center pencil drawer. All the drawers lock, except the pencil drawer. The knee space provided is 26"W.



Item #	Size
1131K	60"W x 30"D x 30"H

instructor's desks



5' INSTRUCTOR'S DESK—Complete demonstration center constructed of solid oak, oak hardwood and oak veneers. It is finished with a chemical resistant, earth-friendly UV finish. Included are four drawers (three small, one file), and a large cupboard. The drawers are expertly crafted from the finest materials with dovetail joinery throughout. Drawers run on epoxy-coated steel glides with nylon bearings for smooth operation. All units come with a GFI protected AC duplex receptacle. You can order with/without these options: combination water/gas fixture, and polypropylene trap. Comes standard with pull-out writing board, support rod with crossbar, rubber base molding and locks for cupboard and all drawers. Dimensions;

- Knee space: 22"W x 32-1/2"D
- Pull-out: 3/4"H x 16-1/4"W x 25-1/4"D
- Overall: 60"L x 30"D x 36"H.
- Small Drawers: 13-1/2"W x 18"D x 4"H
- File Drawer: 13-1/2"W x 18"D x 10-1/4"H
- Cupboard: 19-1/2"W x 28"D x 29-1/2"H
- Sink: 14"W x 10"D x 6"H



1216K-L

Item #	Sink & Fixtures	Top	Sink Side
1214K-L	Yes	Solid Phenolic	Left
1214K-R	Yes	Solid Phenolic	Right
1216K-L*	Yes	Solid Epoxy	Left
1216K-R	Yes	Solid Epoxy	Right
1214KF-L	No	Solid Phenolic	Left
1214KF-R	No	Solid Phenolic	Right
1216KF-L	No	Solid Epoxy	Left
1216KF-R	No	Solid Epoxy	Right



ADA INSTRUCTOR'S DESK

—An accessible instructor's desk. This unit has a roll-under sink with paddle handle faucet and gas fixtures. The knee space is wide enough for a wheel chair to fit into the space. The 34"H meets ADA requirements (having a 19" depth, 27" knee height, 30" width and 8-1/2" knee clearance from the front table top edge). It includes a duplex outlet, file drawer, three smaller drawers and a pull-out writing board. It is finished with a chemical resistant, earth-friendly UV finish. Choose from either a 1" solid epoxy or phenolic resin top with an epoxy sink and a rod socket. Unit also includes a greenlaw arm. Overall size: 72"W x 30"D x 34"H.

Item #	Top
1214K-L-ADA	Solid Phenolic
1216K-L-ADA	Solid Epoxy



accessories




SIDE DESK FOR 5' AND 8' INSTRUCTOR'S DESK—It attaches to the Instructor's Desk. Crafted from solid oak, oak and hardwood veneers (same as the instructor's desk). It is finished with a chemical resistant, earth-friendly UV finish. The side desk features one 13-1/2"W x 18"D x 10-1/4"H file drawer, two 13-1/2"W x 18"D x 4"H storage drawers and one 19-1/2"W x 18"D x 3"H pencil drawer. All drawers come with locks. The unit includes a 14-1/2"W x 20"D pull-out writing board. The classroom side of the unit features a bookshelf with one fixed shelf. The unit can be attached to either the right or left side of the Instructor's Desk.

Item #	Size	Top
1144K	48"W x 30"D x 30"H	Solid Phenolic
1146K	48"W x 30"D x 30"H	Solid Epoxy

LABHAND—Designed to provide a safe, convenient surface for laptop computers in a laboratory setting. The LabHand is mounted on a solid plastic rod which fits securely into standard rod sockets. The platform is made of silver colored steel to create a surface that will withstand the typical chemicals found in laboratories. The LabHand holds laptops and books off the work surface protecting them from spills and increasing the amount of workspace available. 12" square with 5-1/2" rod.

Item #	Product
100994*	LabHand
100994X6*	LabHand (Set of 6)

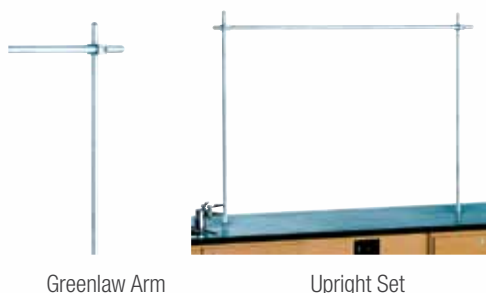

Greenlaw Arm

Upright Set



GREENLAW ARM—Includes 3/4" x 3/4" clamp, a 3/4"x36" upright and a 3/4" x 18" crossbar. Fits in a rod socket with a 3/4" receiver.

UPRIGHT SET—Includes two each 3/4" x 1/2" clamp, two 3/4" x 36" upright and a 1/2" x 40" crossbar. Fits in a rod socket with a 3/4" receiver.

Item #	Product
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar

UPRIGHTS, CROSSBARS, CLAMPS—Solid machined aluminum, these units are nearly impossible to break.

Item #	Size	Product
100002*	3/4" x 3/4"	Clamp
100003*	3/4" x 1/2"	Clamp
100004*	3/4" x 36"	Upright

Item #	Size	Product
100005*	3/4" x 18"	Crossbar
100006*	3/4" x 48"	Crossbar
100007*	1/2" x 40"	Crossbar



accessories



Item #	Product
206505 *	Rod Socket



ROD SOCKETS—Receives the tapered end of 3/4" uprights.



Item #	Product
100047*	Single Gas Cock
100048*	90° Gas Cock
100049*	180° Gas Cock
100466*	360° Gas Cock



GAS COCKS—These fittings have a ground key valve for long life. Units are complete with nipple, lock nut and washer for installation and are deck mounted.



Item #	Product
100074*	Multiservice Fixture
100051*	Cold Water Faucet
100061*	Hot Water Faucet
100467*	Mixing Faucet



GOOSENECK FAUCETS—Chrome plated brass, deck mounted faucets are 9-1/2" H. Multiservice fixtures includes water faucet and two gas cocks. Mounting hardware provided.



Item #	Product
100013*	Chrome Hand Pump



HAND PUMP—Piston action mechanism for high water flow and long term prime lock. Chrome base and spout.



Item #	Product
100728X4*	Floor Mounting Brackets (Set of 4)



FLOOR MOUNTING BRACKETS—Set of four L-shaped mounting brackets are used to permanently attach tables to the floor.

Item #	Product
100118*	2-1/2 Gallon Jug



2-1/2 GALLON JUG—Plastic jug found in our mobile sink units. Holds up to 2-1/2 gallons of liquid.



Item #	Product
100076X4*	Leg Boots (Set of 4)



RUBBER LEG BOOTS—Fit snugly to 2-1/4" table legs. Sold in set of 4.





mobile demo units



Fold-down Monitor Display



Monitor w/Plexiglass



Protect Monitor
w/Phenolic Panel

new!

INTERACTIVE MOBILE MEDIA CABINET—A unique mobile cabinet that allows the user to store, move and use many electronic devices at one time. The unit has ample storage with two dovetailed drawers and cabinets with adjustable shelves for things like projectors, printers, laptops, PCUs and general supply items. The 1" phenolic resin top is both chemical resistant and easy to clean. The top has a cutout section (21.75"W x 17.25"D) where a monitor can be attached to a flat panel monitor arm. A monitor (not included) can adjust to be exposed above the cabinet, or used in the cabinet and seen through a Plexiglas panel that fits over the opening. A phenolic panel is also included for times when the monitor is not needed. The unit includes a pullout keyboard tray, 2 data ports, an extension cord, a duplex outlet and an inlet plug. On the backside of the unit, doors provide access to a power-strip, which has 8 outlets. Throughout the unit there are holes for running cables to the appropriate areas. Venting is cleverly disguised in the bottom of the unit. This unit is available in either oak or maple veneer and is finished in a chemical resistant finish. Overall size: 72"W x 28"D x 36"H.



Item #	Product
4544K	Oak Mobile Media Cabinet
4544M	Maple Mobile Media Cabinet

Note: Monitor, computer, etc. are not included.

mobile demo units



EXTRA LARGE MOBILE DEMO CENTER—Self-contained demonstration center is constructed of first quality solid hardwood, and oak veneers with a chemical resistant, earth-friendly UV finish. This unit features one file drawer with full extension glides, 3 smaller drawers (all dovetailed) and two cupboard spaces all with locks and wire pulls. The drawers run smoothly on epoxy-coated steel glides with nylon rollers. This demonstration center rolls easily on 4" locking swivel casters. Included are 15" x 13"x 5" stainless steel sink, hand water pump, two 2-1/2 gallon water bottles and connecting hose, two rod sockets, two each 3/4" diameter x 36"L upright anodized aluminum rods, and 1 each 1/2" diameter anodized aluminum crossbar with connecting clamps. Included to help organize storage is a handy removable storage tote tray and a pegboard with hooks. Also included is the oak framed combination mirror/marker board which features a shatter-proof plastic mirror and can be reversed to serve as a handy marker board. The combination mirror/marker board can be conveniently stored in one of the cupboard spaces. Electrical service is supplied by means of a 25' extension cord which can be connected to a GFI protected AC duplex receptacle installed in the side of the cabinet. The top is black 1-1/4" ChemGuard laminate. Unit can also be ordered without the sink, faucet, rod sockets, upright set and mirror/marker board. Mobile capacity: 500 lbs. Overall dimensions: 30"D x 54"W x 36"H. Shipping weight 290 lbs.



Item #	Product
4342K*	Mobile Demo Unit w/Sink & Fixtures
4342KF	Mobile Demo Unit w/Flat Top
4342KF-RS	Mobile Demo Unit w/Rod Sockets
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar

NOTE: ON "RS" UNITS, THE GREENLAW ARM OR UPRIGHT FULL CROSSBAR SET IS NOT INCLUDED.



MOBILE DEMO CART—Designed for the transportation of equipment and supplies. Features solid oak construction. The legs are 2-1/4" square and are securely fastened to the 4-1/2" high apron with factory installed 5/16" hanger bolts through the heavy-duty 13 gauge steel corner braces. Includes leg stretchers and auxiliary oak plywood shelf. The unit rolls easily on 4 factory installed 3" ball-bearing swivel casters with brakes. The cart is finished with a chemical resistant, earth-friendly UV finish. The unit is supplied with either a factory installed black 1-1/4" Plastic Laminate or 1-1/4" chemical resistant ChemGuard top. Shipped partially unassembled. Dimensions: 36"W x 24"D x 30"H. Wt. 95 lbs.

Item #	Top
4501K*	P. Lam.
4502K*	ChemGuard



mobile demo units



USB Serial Port
Connection/S-Video/RCA Plug

MOBILE DESK WITH CAMERA MOUNT—Designed to hold and lock the Ken-A-Vision Furniture Cam. The mobile desk is equipped with two storage compartments to hold your camera, and other large items. The three drawers (one file, two box) have dovetail joints and run smoothly on epoxy-coated steel glides. Includes pull-out writing board and all drawers and doors have locks. This unit features a 1-1/4" chemical resistant ChemGuard top, 21" x 28" pegboard and assorted hooks for vertical storage, 15" x 13" x 5" stainless steel sink and hand water pump. Also comes with two 2-1/2 gallon water bottles with all connecting hoses, two 3/4" diameter upright rods, one 1/2" diameter crossbar, two connecting clamps and two rod sockets installed in the table top. A GFI protected AC duplex electrical outlet is built into the lab, with an additional outlet for monitor and camera plug-in. Unit has four heavy-duty 4" locking swivel casters. There is a storage area for the optional mirror/marker board. It is finished with a chemical resistant, earth-friendly UV finish. Dimensions: 28"D x 48"W x 36"H. Weight: 235 lbs.

Item #	Product
4332KAV	Mobile Desk w/Camera Mount
4000K*	Optional Mirror
4001K*	Optional Mirror/Marker Board

MONITOR AND CAMERA MOUNT NOT INCLUDED

MOBILE DESK WITH STORAGE—Same unit as above without the camera mount. Mobile capacity: 500 lbs. Unit can also be ordered without the sink, faucet, rod sockets and upright set. Weight: 235 lbs.



Item #	Product
4332K*	Mobile Desk w/Storage
4332KF	Mobile Desk w/Storage and Flat Top
4332KF-RS	Mobile Desk w/Storage Flat Top and Rod Sockets
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar
4000K*	Optional Mirror
4001K*	Optional Mirror/Marker Board

NOTE: ON "RS" UNITS, THE GREENLAW ARM OR UPRIGHT FULL CROSSBAR SET IS NOT INCLUDED.



4332K

mobile demo units



MOBILE LABORATORY UNIT—A fully equipped, locking self-contained demonstration unit. It features, a sturdy 1-1/4" thick high-pressure laminated surface with durable drop leaf extensions. Constructed from select oak and oak veneers with a chemical resistant, earth-friendly UV finish. It features unbreakable aluminum sliding door glides and lock. The unit comes standard with a stainless steel sink, and hand pump. A GFI protected AC duplex electrical outlet is built into the lab. The cabinet rides on 4" heavy-duty ball-bearing swivel casters with brakes. Includes a 25- foot extension cord, two tote trays, two 2-1/2 gallon water bottles, connecting hoses, two 3/4" diameter upright rods, one 1/2" diameter crossbar and clamps, a propane burner with stand and 2 adjustable shelves. Unit is also available without the sink, faucet, rod sockets and upright set. Mobile capacity: 500 lbs. Weight: 230 lbs. Overall dimensions: 48"W x 24"D x 36"H. Extended Top: 72"W x 24"D. Sink: 15"W x 13"L x 5"D.



4121K

Item #	Product
4121K *	Mobile Laboratory Unit
4121KF	Mobile Laboratory Unit w/ Flat Top
4121KF-RS	Mobile Laboratory Unit Mobile w/Flat Top, and Rod Sockets
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar
4000K*	Optional Mirror
4001K*	Optional Mirror/Marker Board

NOTE: ON "RS" UNITS, THE GREENLAW ARM OR UPRIGHT FULL CROSSBAR SET IS NOT INCLUDED.



4222K



MOBILE DEMO TABLE WITH DRAWERS—A self-contained unit with its own water supply, waste container, stainless steel sink and hand pump, two 3/4" diameter upright rods, one 1/2" diameter crossbar and clamps. The four drawers are mounted on epoxy-coated nylon bearing glides. Sliding, locking doors conceal the drawers and storage compartment, and the supplied two 2-1/2 gallon water bottles. It is finished with a chemical resistant, earth-friendly UV finish. Heavy-duty 4" ball-bearing swivel casters with brakes, support 500 lbs. Includes a GFI protected AC duplex electrical outlet and can be powered with the supplied 25' extension cord. The top is 1-1/4" thick ChemGuard high pressure laminate. Sink measures 15"W x 13"L x 5"D. Unit is also available without the sink, faucet, rod sockets and upright set. Overall dimensions are 48"W x 28"D x 36"H. Weight: 215 lbs.

Item #	Product
4222K*	Mobile Demo Table
4222KF	Mobile Demo Table w/Flat Top
4222KF-RS	Mobile Demo Table w/Flat Top and Rod Sockets
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar
4000K*	Optional Mirror
4001K*	Optional Mirror/Marker Board



mobile demo units



4111K



ECONOMY MOBILE LAB TABLE—Constructed from select oak and oak veneers with a chemical resistant, earth-friendly UV finish. It features unbreakable aluminum sliding door glides and lock. Edges of the 1-1/4" thick high-pressure laminated top are bound with a Plastic-T Molding. With a mobile capacity of 500 lbs., it is equipped with a stainless steel sink, hand pump, heavy-duty 4" locking swivel casters, 25' grounded extension cord, two 2-1/2 gallon water bottles, connecting hoses, two 3/4" diameter upright rods, one 1/2" diameter crossbar with clamps, two fully adjustable shelves, and a GFI protected AC duplex electrical outlet. Unit is also available without the sink, faucet, rod sockets and upright set. Sink measures 15"W x 13"L x 5"D. Dimensions: 48"W x 24"D x 36"H. Weight: 185 lbs.

Item #	Product
4111K*	Economy Lab Table
4111KF	Economy Lab Table w/Flat Top
4111KF-RS	Economy Lab Table w/Flat Top, and Rod Sockets
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar
4000K*	Optional Mirror
4001K*	Optional Mirror/Marker Board

NOTE: ON "RS" UNITS, THE GREENLAW ARM OR UPRIGHT FULL CROSSBAR SET IS NOT INCLUDED.

ADA MOBILE LAB UNIT—Designed as a work station for the physically impaired or for classroom demos. Constructed of solid oak, hardwood and oak veneers with a chemical resistant, earth-friendly UV finish. The 1-1/4" ChemGuard top is resistant to most lab chemicals. It is supplied with or without a stainless steel sink, toggle switch water pump, two 2-1/2 gallon supply and waste containers, GFI AC duplex receptacle with inlet socket, 25' extension cord, two rod sockets, two each 3/4" diameter x 36" long upright rods, 1/2" x 40" crossbar and two rod clamps. With a mobile capacity of 500 lbs. it rolls easily on 4" ball-bearing swivel casters with brakes. Dovetailed drawers are constructed of plywood and ride on epoxy coated glides. Cupboard door and drawers have locks. Can be accessed with a wheelchair from two sides. May be used with optional demonstration mirror. Sink measures: 15"L x 13"W x 5"D. Inside drawers measure: 14.5"W x 19"D x 6"H. Dimensions: 72"L x 24"D x 33"H. Wheelchair access: 32"W x 29"H. Weight: 240 lbs.

Item #	Product
4352K*	ADA Compatible Mobile Lab Station
4352KF	ADA Compatible Mobile Lab Station w/Flat Top
4352KF-RS	ADA Compatible Mobile Lab Station w/Rod Sockets
100245GRARM*	Greenlaw Arm
100604*	Upright Set w/Full Crossbar
4000K*	Optional Mirror
4001K*	Optional Mirror/Marker Board

NOTE: ON "RS" UNITS, THE GREENLAW ARM OR UPRIGHT FULL CROSSBAR SET IS NOT INCLUDED.



4352K



mobile storage

MIRROR—The shatter resistant acrylic mirror is enclosed in a solid oak frame with a natural transparent finish. Easily mounts on 1/2" diameter crossbar of any of our mobile labs which include uprights and crossbar. It can be tilted to any angle—held secure with nylon screws in the mounting brackets. Measures 22-1/2" x 34-1/2". Weight 10 lbs.

Item #	Product
4000K*	Mirror



MIRROR/MARKER BOARD—This solid oak frame unit allows for the unobstructed use of the shatter resistant acrylic mirror and marker board. Unit can be tilted to any angle or completely reversed. The unit is held secure to the 1/2" diameter crossbar (not included, see page 22) by means of screws in the frame. Measures 20-3/4" x 27-3/4". Weight 12 lbs.

Item #	Product
4001K*	Mirror/Marker Board



4101K



BARE BONES MOBILE STORAGE CABINET—Ideal for moving demonstrations around the room, school or lab. Constructed of select oak and oak veneers with a chemical resistant, earth-friendly UV finish. Includes unbreakable aluminum sliding door glides, lock and one fixed shelf. Edges of the 1-1/4" thick high pressure laminated top are bound with a Plastic-T Molding. With a weight capacity of 500 lbs., it rolls on 4" heavy-duty ball-bearing swivel casters with brakes. Dimensions: 48"W x 24"D x 36"H. Weight: 180 lbs.

Item #	Product
4101K*	Bare Bones Mobile Storage Cabinet

MOBILE STORAGE CABINET—A compact, versatile rolling storage cabinet that is constructed of select oak and oak veneers with a chemical resistant, earth-friendly UV finish. The adjustable shelf is arranged to allow tall items to be placed in the unit as well as allowing use of the shelf. With a mobile capacity of 500 lbs., it comes equipped with four heavy-duty swivel casters with brakes. The doors lock for security. The unit is available with black 1-1/4" Plastic Laminate or 1-1/4" ChemGuard top. Dimensions: 36"W x 24"D x 36"H. Shipping weight: 140 lbs.

Item #	Top
4401K*	P. Lam
4402K*	ChemGuard



4401K





mobile storage



MOBILE TOTE TRAY STORAGE—Store various materials with these large totes. Constructed of oak veneers and solid oak with a chemical resistant, earth-friendly UV finish. With a mobile capacity of 500 lbs., this unit rolls easily on 4" locking swivel casters. The two lockable doors swing 270° and conceal 21 tote trays. Tote trays measure 19"L x 14-1/4"W x 3-1/2"H. Plastic laminate top can double as a work surface. Overall unit measures: 48"W x 24"D x 41-1/2"H.

Item #	Product
4751K*	Mobile Tote Tray Storage
100135*	Replacement Tote Tray

WRITE-N-ROLL MOBILE STORAGE—Versatile mobile storage with an adjustable position top that has a 3mm black PVC edge. The top is made of 3/4" white marker board so messages or techniques can be demonstrated. Directly beneath the top is storage for flat items like paper and books measuring 19" x 25." The top can be positioned conveniently with the two 10-position arms. The cabinet is constructed of oak and select hardwood veneers with a chemical resistant, earth-friendly UV finish. Open the doors which feature full-length aluminum extrusion handles to reveal even more storage. The left side of the cabinet has five tote trays (13" x 17-1/2" x 3-1/2" inside dimension) that can be adjusted on 1-1/2" centers. The right side of the unit features three drawers (13" x 17" x 2-1/2" inside dimension) that are mounted on adjustable drawer glides that can be mounted on 1/2" centers. With a mobile capacity of 500 lbs., it rolls easily on 4" swivel casters with brakes. Dimensions: 36"W x 24"D x 36"H. Weight: 160 lbs.



Item #	Product
4901K	Write-N-Roll Mobile Storage



MOBILE BALANCE STORAGE CABINET—Constructed of solid oak and oak veneered plywood with chemical resistant, earth-friendly UV finish. Featured are locking doors, both front and back with unbreakable aluminum glides. Supplied with four adjustable shelves (weight capacity of 40 lbs per sq. ft.). With a mobile capacity of 500 lbs., it rolls on heavy-duty swivel casters with brakes. The top is black 1-1/4" plastic laminate and furnishes an auxiliary work surface. Dimensions; Shelves are 3/4"H x 8"D x 42"L, Overall is 48"W x 24"D x 40"H. Weight: 230 lbs.

Item #	Product
5201K*	Mobile Balance Cabinet Storage

mobile storage



MOBILE MICRO-CHARGE STATION—A mobile oak cabinet that is equipped to recharge any **Ken-A-Vision** cordless microscope within eight hours while it is safely stored behind locked doors. The peninsula style cabinet interior (7"W x 9"D x 15"H) holds 16 microscopes and has storage space for other essentials such as slides. The LED meter lets you know when the microscopes are charging and when they are fully charged. All this with only one electrical connection from the cabinet's side. Is is finished with a chemical resistant, earth-friendly UV finish. The 1-1/4" plastic laminate top doubles as a work surface. Mobile capacity: 500 lbs. Dimensions: 48"W x 30"D x 41-1/2"H. Weight: 311 lbs.

Item #	Product
4741K*	Mobile Micro-Charge Station



MOBILE MICROSCOPE STORAGE CABINET—Provides a safe, convenient way to store your microscopes. It is constructed of solid hardwood, solid oak and oak veneers with a chemical resistant, earth-friendly UV finish. The front and back doors have locks. The unit contains removable wire dividers to provide an individual space for each individual microscope (6" x 11" x 16"). With a mobile capacity of 500 lbs., this cabinet rolls on heavy-duty swivel casters with brakes. The top is 1-1/4" black plastic laminate and furnishes an auxiliary work surface. The cabinet has storage for 24 standard classroom microscopes. Dimensions: 48"W x 24"D x 40"H. Weight: 202 lbs.

Item #	Product
4701K*	Mobile Microscope Storage Cabinet




MOBILE LAPTOP RECHARGING STATION—Sliding hardwood plywood doors front and back provide safe, easy access to the charging station. Multiple outlet plug strips in the back section allow easy plug-in of up to 20 computers. The cords from the outlet strips can be fed through a hole in the floor of the unit so that the doors can be locked during charging. The wire dividers minimize heat build up during charging. The unit is also equipped with vents and fan to disperse heat. The unit is constructed of oak, oak plywood and hardwood plywood with a 1-1/4" black plastic laminate top. With a mobile capacity of 500 lbs., this cabinet rolls on 4" heavy duty swivel locking casters. Laptop opening: 3-5/8"W x 13"D x 13"H. Dimensions: 48"W x 24"D x 36"H.

Item #	Product
5301K	Mobile Laptop Storage Station





tall storage

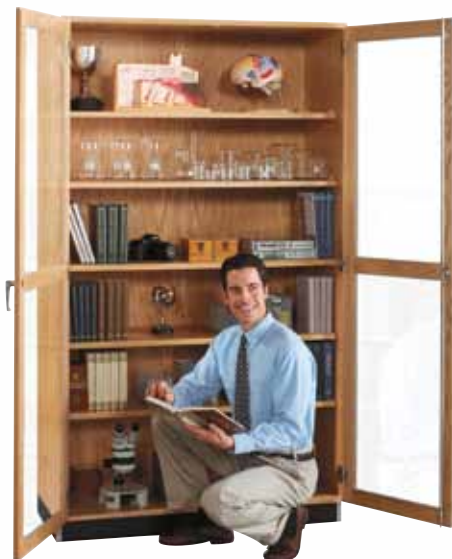

TALL STORAGE W/GLASS DOORS—Showcase your specimens with this beautiful clear glass door cabinet. Solid hardwood, solid oak and oak veneers, protected by a chemical resistant, earth-friendly UV finish. The oak framed 3/16" tempered glass doors are designed to crumble when broken preventing large dangerous shards. Each cabinet comes with a three-point locking system. The 24"W cabinet has one locking door that is hinged right. Glass door units have two fixed and four adjustable 1" shelves (weight capacity: 40 lbs. per sq. ft). Dimensions: 22"D x 84"H.



Item #	Product
315-2422	24"W Glass Door Tall Storage
358-3622	36"W Glass Door Tall Storage
358-4822*	48"W Glass Door Tall Storage

Cabinet should be secured to the wall for safety. Mounting hardware is not included.

TALL STORAGE W/OAK VENEER DOORS—Solid hardwood, solid oak and oak veneers with a chemical resistant, earth-friendly UV finish. Available with hinged oak doors. Each cabinet comes with a three-point locking system. Units 36" and wider have 1" shelves. 24" wide oak door units have 3/4" shelves, and one locking door that is hinged right. All units include two fixed and four adjustable shelves (weight capacity: 40 lbs per sq. ft.)



Item #	Product
313-2422	24"W Oak Door Tall Storage
353-3622	36"W Oak Door Tall Storage
353-4822*	48"W Oak Door Tall Storage

Cabinet should be secured to the wall for safety. Mounting hardware is not included.



TOTE TRAY STORAGE CABINET—This organized storage case holds 48 tote trays (tray dimensions: 19"W x 14-1/2"D x 3-1/2"H). A great place for students to store their art or science projects. Tote tray storage case has hinged oak doors with three-point locking handles. It is finished with a chemical resistant, earth-friendly UV finish. Weight: 450 lbs. Overall Dimensions: 48"W x 22"D x 84"H.



Item #	Product
350-4822*	Tote Tray Storage Cabinet
100135*	Replacement Tote Tray

Cabinet should be secured to the wall for safety. Mounting hardware is not included.

tall storage



TOTE TRAY & SHELVING STORAGE CABINET—A great place to store projects and larger objects...perfect for a science or art room. This unit is made up of solid oak and oak veneers with an earth-friendly, chemical resistant UV finish. The doors have three-point locking handles for secure storage. Inside you will find 24 totes (19"W x 14-1/2"D x 3-1/2"H) and one fixed and two adjustable shelves. Dimensions: 48"W x 22"D x 84"H.

Item #	Product
351-4822	Tote Tray & Shelving Storage Cabinet
100135*	Replacement Tote Tray

Cabinet should be secured to the wall for safety. Mounting hardware is not included.



ROCK/PAPER STORAGE CABINET—Whether you are storing paper, rocks, or insects this unit will fit the bill. Seven large drawers (4"H x 44"W x 26"D) big enough to hold large flat objects. The three shelves (two adjustable) allow you to store over-sized items like models. Unit is constructed of oak and oak veneers with an earth-friendly, chemical-resistant UV finish. It has a three-point locking handle system. Overall Dimensions: 48"W x 30"D x 84"H.

Item #	Product
354-4830	Rock/Paper Storage Cabinet

Cabinet should be secured to the wall for safety. Mounting hardware is not included.



MICROSCOPE STORAGE CABINET—Constructed of solid oak, hardwood and oak veneers, using dowel construction, with a chemical resistant, earth-friendly UV finish. Solid oak framed glazed doors are hinged with three-point locking handles. Tempered glass is 3/16" thick. Microscopes are stored safely and conveniently on two fixed and three adjustable 1" shelves (weight capacity of 40 lbs per sq. ft.). Adjustable wire dividers allow for storage of microscopes of various sizes. Dividers can be removed to use shelves for general storage. The 24"W cabinet has one locking door that is hinged right. Available in three sizes. 24"W x 16"D x 84"H holds 15 microscopes. 36"W x 16"D x 84"H holds 20 microscopes. 48"W x 16"D x 84"H holds 30 microscopes. Individual compartment measures 7-1/4"W x 14"D x 14-1/2"H.

Item #	Size/Weight
370-2416	24"Wide/Wt. 230 lbs.
372-3616	36" Wide/Wt. 255 lbs.
372-4816*	48" Wide/Wt. 326 lbs.





tall storage



MICRO-CHARGER GLASS DOOR CABINET—This beautiful oak cabinet with tempered glass doors is finished with a chemical resistant, earth-friendly UV finish. It is equipped to recharge up to 16 **Ken-A-Vision** cordless microscopes (not included) within eight hours while it is safely stored behind locked doors (Microscope Dimensions: 8-3/4"W x 12-1/4"D x 16-3/4"H). The LED meter lets you know when the microscopes are charging and when they are fully charged. All this with only one electrical connection. The microscopes are protected from theft by a three-point locking handle. Overall size: 36"W x 16"D x 84"H.

Item #	Size
373-3616	36"W x 16"D x 84"H

Cabinet should be secured to the wall for safety. Mounting hardware is not included.

WARDROBE CABINET—Constructed of oak veneers and solid hardwoods with a chemical resistant, earth-friendly UV finish. Half of the unit features a coat rod with a shelf above. Unit also includes two unexposed file drawers with full extension glides. One fixed and three adjustable shelves are found above the file drawers. Mirror measures: 13"H x 9.5"W. Unit has three-point latching handles which lock. Overall dimensions: 36"W x 22"D x 84"H. Weight: 360 lbs.

Item #	Product
360-3622	Wardrobe Unit

Cabinet should be secured to the wall for safety. Mounting hardware is not included.



SKELETON CABINET—Choose the model that suits your classroom skeleton—either on a rolling stand or a hanging model on a retractable glide. Constructed of solid oak and oak veneers, with a chemical resistant, earth-friendly UV finish. The retractable glide is included in the hanging model. Both cabinets have wood framed 3/16" tempered glass doors for easy viewing. The cabinets can be locked for security purposes. A shelf is included with each unit for extra storage. Overall size: 24"W x 22"D x 84"H. Wt. 220 lbs.

Item #	Product
375-2422	Hanging Model Cabinet
377-2422	Roll-in Cabinet

*Cabinet should be secured to the wall for safety. Mounting hardware is not included.
SKELETON NOT INCLUDED*

shower/eyewash stations



SAFETY STATION—Complete all-in-one safety station has a built-in shower and eyewash station, and storage for other equipment like fire blankets, spill kits, first aid kits and a fire extinguisher. The shower head has a high volume continuous water flow. Pulling the conveniently located pull handle activates the stay open valve. The eyewash is mounted on black plastic laminate over a stainless steel bowl but can be detached and pulled up to 12' away from the station for washing eyes or other small areas. The eyewash paddle handle activates the stay open valve leaving the users hands free. Shelving, both locking and open, is available. The area beneath the storage is open to allow the mounting of various safety apparatus. The unit is ADA compliant with room for a wheelchair and lever controlled fixtures. It is finished with a chemical resistant, earth-friendly UV finish. Dimensions: 48"W x 22"D x 84"H. Weight: 280 lbs.

Item #	Product
6950K	Safety Station



SHOWER/EYEWASH STATION—This unit provides convenient storage for lab safety supplies. Featured is an emergency drench shower head of high impact ABS plastic with a stay open chrome-plated brass valve. Eyewash unit has a 10" diameter stainless steel bowl, lever action ball valve on a plastic laminate top. Upper storage space has 2 adjustable shelves. The lower cabinet locks for security or access for a wheelchair on ADA model. Service piping not supplied. It is finished with a chemical resistant, earth-friendly UV finish. Dimensions: 36"W x 22"D x 84"H. Weight: 265 lbs.

Item #	Product
6900K	Shower/Eye Wash
6900K-ADA	Shower/Eye Wash ADA Model





casework



121-3622 | Base Cabinet



106-3622 | Base Cabinet



103-3622 | Base Cabinet

BASE CABINETS—Constructed of solid hardwood, solid oak and oak veneers. Base cabinets with choice of three styles: door cabinet, door/drawer cabinet or drawer cabinet. Cabinets can be ganged together to create a wall of base cabinets. All cabinets measure 36"W x 22"D x 35"H. They are finished with a chemical resistant, earth-friendly UV finish. Tops and rubber base molding sold separately. Base molding sold by the foot.

SOLID RESIN TOPS—Finish off your base cabinets with these solid epoxy resin or phenolic resin tops. These solid tops provide superior chemical resistance. Tops include 4" back splash. Available in 1" thick black.

WALL CABINETS—Constructed of solid hardwood, solid oak and oak plywood veneers. Choose from two models: Oak veneered solid doors or solid oak framed glass doors. Glass doors feature 3/16" thick safety glass. Cabinets have a lock and two adjustable shelves. They are finished with a chemical resistant, earth-friendly UV finish. Units come in varying widths and are 12"D x 30"H.



D06-3612 | Wall Cabinet



D03-3612 | Wall Cabinet

Item #	Product
121-3622	Drawer Base Cabinet
106-3622	Door/Drawer Base Cabinet
103-3622	Door Base Cabinet
D03-3012	30"W Solid Double Doors - Wall
D03-3612	36"W Solid Double Doors - Wall
D03-4212	42"W Solid Double Doors - Wall
D03-4812	48"W Solid Double Doors - Wall
D06-3012	30"W Glass Double Doors - Wall

Item #	Product
D06-3612	36"W Glass Double Doors - Wall
D06-4212	42"W Glass Double Doors - Wall
D06-4812	48"W Glass Double Doors - Wall
207529F	38"W x 24"D Solid Epoxy Top
207530F	74"W x 24"D Solid Epoxy Top
214515F	38"W x 24"D Solid Phenolic Top
214516F	74"W x 24"D Solid Phenolic Top
100283*	Rubber Base Molding per foot

book trucks



BT012



BT1113



BT311



Northwoods Oak



Rich Honey Oak



Deep Bronze Oak

SOLID OAK BOOK TRUCKS—Bring the beauty of wood into your classroom or lab. Durable trucks combine a distinctive wood grain of select oak with a natural, medium or dark finish. The fixed shelves are glued then reinforced with screws to support heavy volumes. All trucks have 5" locking swivel casters.

Item #	Color	# of Shelves	Clearance	Size	Wt.
BT011	Northwoods Oak	2 Sloped and 1 Flat	12"	37-3/4" x 32" x 13"	60 lbs.
BT012	Rich Honey Oak	2 Sloped and 1 Flat	12"	37-3/4" x 32" x 13"	60 lbs.
BT013	Deep Bronze Oak	2 Sloped and 1 Flat	12"	37-3/4" x 32" x 13"	60 lbs.
BT111	Northwoods Oak	3 Flat	12"	37-3/4" x 32" x 13"	55 lbs.
BT112	Rich Honey Oak	3 Flat	12"	37-3/4" x 32" x 13"	55 lbs.
BT113	Deep Bronze Oak	3 Flat	12"	37-3/4" x 32" x 13"	55 lbs.
BT211	Northwoods Oak	3 Flat	15"	45" x 40-1/2" x 14"	70 lbs.
BT212	Rich Honey Oak	3 Flat	15"	45" x 40-1/2" x 14"	70 lbs.
BT213	Deep Bronze Oak	3 Flat	15"	45" x 40-1/2" x 14"	70 lbs.
BT311	Northwoods Oak	4 Sloped and 1 Flat	11"	42-1/4" x 42" x 18"	162 lbs.
BT312	Rich Honey Oak	4 Sloped and 1 Flat	11"	42-1/4" x 42" x 18"	162 lbs.
BT313	Deep Bronze Oak	4 Sloped and 1 Flat	11"	42-1/4" x 42" x 18"	162 lbs.



bookcases / stools



STORAGE BOOKCASE—Open oak shelving allows for a clear view (two fixed/four adjustable shelves). A chemical resistant, earth-friendly UV finish coats the top, back and sides of each unit. Available in three widths. (

CHEMICAL BOOKCASE—Open oak shelving with added safety, anti-roll lips are features on each fixed shelf (all shelves are fixed). A chemical resistant, earth-friendly UV finish coats the top, back and sides of each unit. Available in three heights.

Item #	Product	Size	Wt.
445-3616	36"H Chemical Bookcase	36"W X 16"D X 36"H	120
446-3616	48"H Chemical Bookcase	36"W X 16"D X 48"H	240
447-3616	72"H Chemical Bookcase	36"W X 16"D X 72"H	310


Item #	Product	Size	Wt.
301-2422	84"H Storage Bookcase	24"W X 22"D X 84"H	250
301-3622	84"H Storage Bookcase	36"W X 22"D X 84"H	270
301-4822	84"H Storage Bookcase	48"W X 22"D X 84"H	310

SOLID OAK STOOLS—Our durable oak stools are available in three sizes that match our lab tables and workstations. They are designed with square post legs, and have a 14" diameter seat.

Item #	Style
5018K	18" Stool
5024K	24" Stool
5030K	30" Stool



cabinet guide

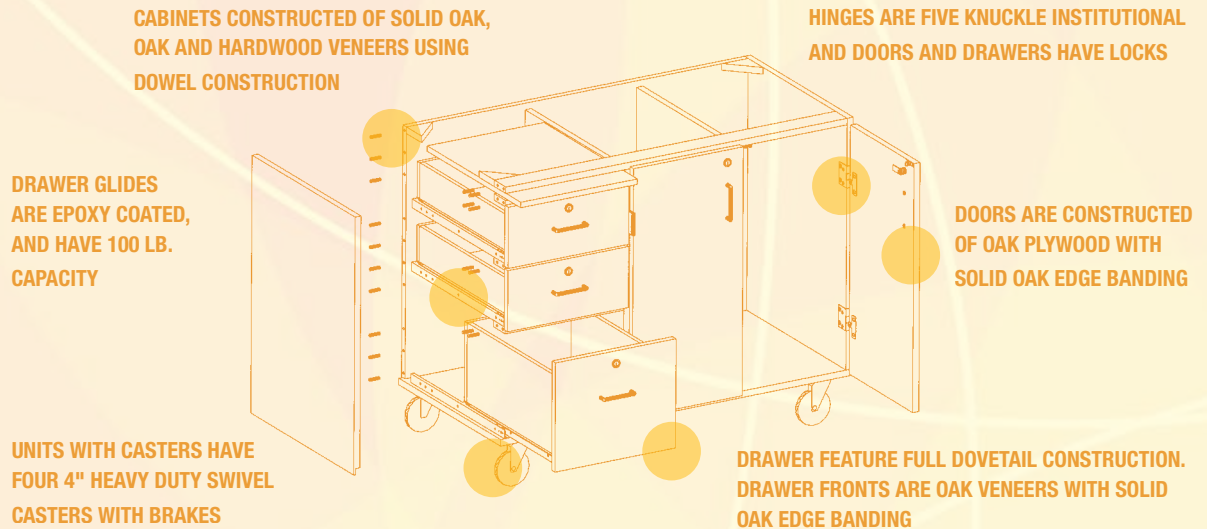

top guide

good

1-1/4" ARP HIGH PRESSURE PLASTIC LAMINATE is widely used for student science tables. It offers an excellent work surface that wears exceptionally well under normal use. Plastic laminate has superior resistance to scratching, with limited resistance to high temperatures and severe chemicals. Top construction consists of permanently bonding a .050 thick sheet of ARP High Pressure Laminate to a 45 lb. density industrial grade particle board core, bottom is covered with a .040" thick brown phenolic backer sheet to properly balance the resulting panel, maximizing it's resistance to warping. Technical data: exceeds 1800 wear cycles and 66 inches impact resistance. Edges of top are finished with vandal resistant T-molding. Corners have 1-3/4" radius. Color: black.

better

1-1/4" CHEMGUARD™ incorporates the decorative features of high pressure laminate into a durable, chemical resistant surfacing material. Tests prove ChemGuard™ has outstanding chemical resistance, plus significantly better NEMA wear value than competitive laminates. ChemGuard™ even has a Class 1 fire rating. A 3mm PVC molding is applied with heat sensitive adhesive to create an attractive waterproof and durable edge. Color: black.

best

1" PHENOLIC RESIN is a composite thermo-fused under heat and pressure to form a solid black chemical resistant composite throughout the entire thickness of the top. It is designed for use as counter tops, pegboards and reagent racks. Phenolic Resin will not stain or etch when exposed to most chemicals or solvents. Phenolic Resin exhibits no effect with subjected to these common laboratory chemicals: Nitric Acid 30%, Acetic Acid 98%, Hydrofluoric Acid 48%, Sulfuric Acid 33%, Methyl Ethyl Ketone, Phosphoric Acid 85% and others. Color: black.

1" SOLID EPOXY RESIN is impervious to normal laboratory chemicals and heat. It's monolithic composition eliminates the possibility of swelling, delamination, and exposure of the substrate. It's soft, low glare surface is simple to clean. Solid Epoxy Resin exhibits no effect when subjected to these common laboratory chemicals: Nitric Acid 70%, Acetic Acid (glacial), Hydrochloric Acid 20%, Sulfuric Acid 60%, Methyl Alcohol, Toluene, Benzene, and others. Solid Epoxy Resin is solid, homogeneous, molded from modified epoxy resins, completely oven cured during processing. It is extremely durable and completely impervious to moisture. Exposed edges and corners are finished with a 1/8" bevel. Color: black.



address

300 S Krueger Street
Suring, WI 54174

phone

877.348.9663
920.842.2136

fax

920.842.2499

e-mail

info@diversifiedwoodcrafts.com

web

www.diversifiedwoodcrafts.com



Mixed Sources

Product group from well-managed
forests, controlled sources and
recycled wood or fiber

Cert no. SW-COC-001886
www.fsc.org
© 1996 Forest Stewardship Council